
Bettgriff "FRIDA"

Art.-Nr. BETTG 010

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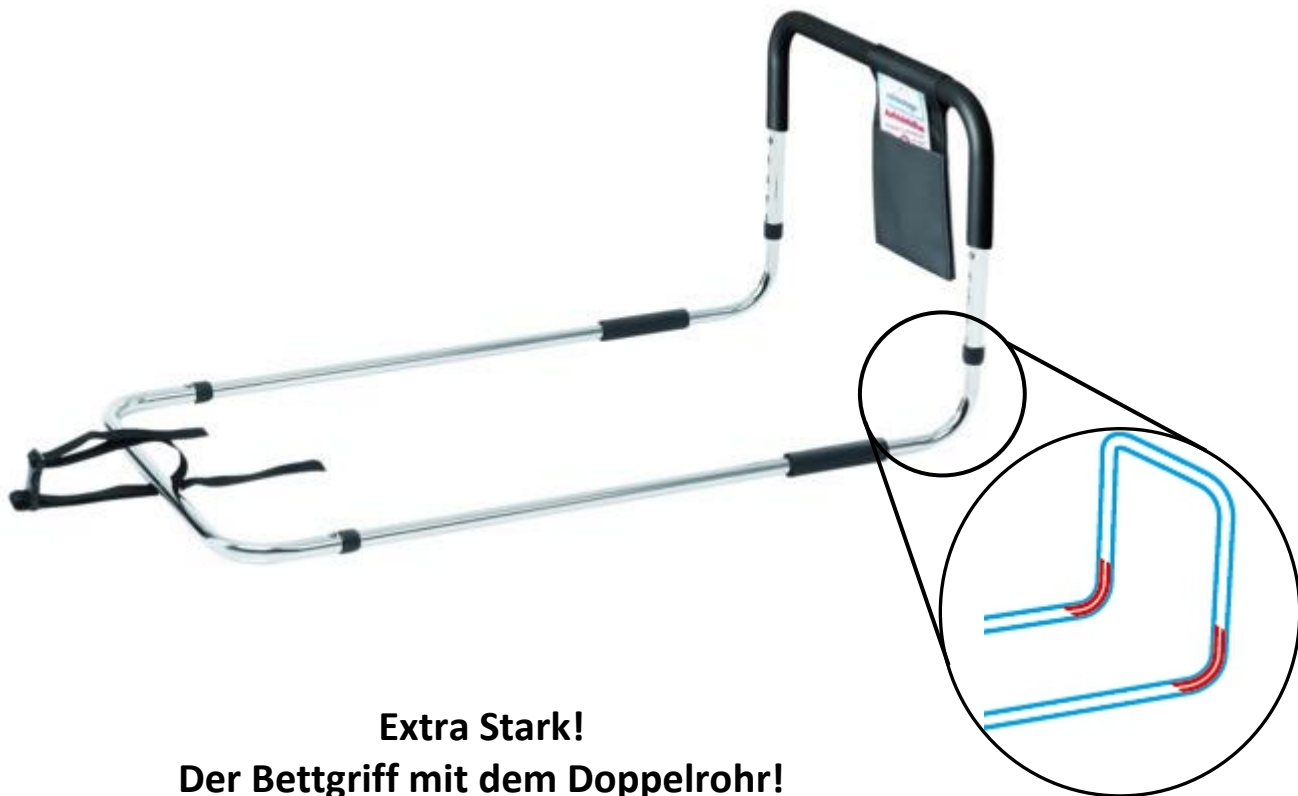1.0 Foreword:

Dear customer,

We are pleased that you have chosen a quality product from *rehastage GmbH*. A good relationship with our customers is of great importance to *rehastage*. Customer proximity means fast service, as little bureaucracy as possible and close cooperation with our customers. Seek advice from your treating physician or therapist to ensure that this product is suitable for you and meets your needs.

⚠ **For your personal safety, please read these instructions carefully before using the product. If there are any problems understanding the information, warnings and instructions listed here, please inquire with your specialist dealer to avoid possible injury or property damage due to incorrect operation.**

Do not throw away these instructions for use even after installing the aid so that you can use it later if necessary.

If you need spare parts or accessories or simply have a question about your aid, we are here for you. We want you to be satisfied with our products and service. We at *rehastage* are therefore constantly working on the further development of our products. For this reason, changes in shape, technology and features may occur in our product range. Therefore, no claims can be derived from the information and illustrations in this instruction manual.

2.0 Intended Use

The use is carried out by the user / patient himself, if necessary, with the support of nursing staff or caring relatives. It allows the user / patient to get into a sitting position without outside help or to support them when getting in and out of bed. The device is intended for use by an adult person/patient who does not suffer from the mentioned contraindications.

- Indication:

The bed handle "FRIDA" (Bed Entry and Exit Aid) is an aid for users and / or patients with generally reduced or impaired movements.

The provision of straightening aids is indicated if, for example, due to a loss of balance, the position between raising the upper body and achieving a safe sitting posture cannot be maintained. For example, at least one arm must be sufficiently loaded and sufficient trunk control and stretchability of the legs must be given to enable the upper body to straighten up largely independently.

- Contraindications:

⚠ The "FRIDA" bed handle is not suitable for people with Alzheimer's or dementia, who take sedatives, are confused or frail. There is an increased risk of entrapment and injury here.

⚠ People who try to climb over the product are at increased risk of injury from falling.

⚠ The device is NOT intended for use by children.

⚠ **Be sure to observe the other warnings and precautions listed under point 5!**

The "FRIDA" bed handle is a Class I medical device in accordance with Regulation (EU) 2017/745 of the European Parliament on medical devices.

3.0 Scope of delivery:

1 piece	Steel tube "Handle part with foam rubber cover"
2 pieces	Steel Pipe "Base Piece Long Side"
1 piece	Steel tube U-shaped, "tail piece"
1 piece	"Belt", 140 cm
1 piece	Bag


4.0 Assembly and assembly instructions for the "FRIDA" bed rail

4.1 Assembling and assembling the bed handle

- 4.1.1 Place the end piece on a flat large surface and insert one of the two base pieces long side into each end until you hear the safety pins snap into place and they also **visibly protrude from the holes on the end piece**. Thus, both base pieces are securely connected to the end piece.

- 4.1.2 Next, put the handle part on the two ends of the base pieces at the same time. Here, too, please check that the safety pins are engaged correctly, which must visibly protrude from the handle part to ensure a secure hold.

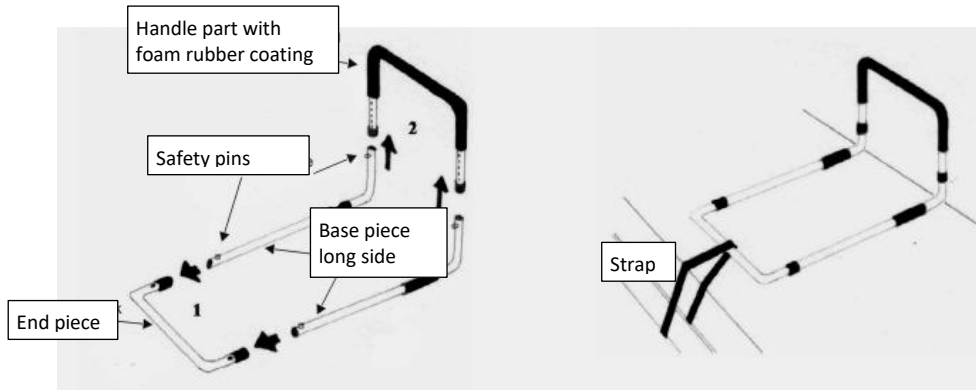4.2 Installation of the bed handle in the bed

- 4.2.1 Lift or remove the mattress and position the "FRIDA" bed handle on the inner bed frame. The handle part should be positioned in such a way that the user can easily reach the handle in a lying position.

⚠ CAUTION: Keep a distance of at least 32 cm between the handle and the head or foot end of the bed to prevent entrapment.
The handle should protrude at least 10 cm above the surface of the mattress to prevent it from climbing over.

- 4.2.2 Place the strap around the end piece of the bed handle and then pull the strap around the opposite frame of the slatted frame. Alternatively, you can also pull the strap around the side part of the bed frame. We recommend attaching it to the slatted frame if possible.

⚠ **Do not use the aid without attaching the seat belt to the bed frame or slatted frame!**

- 4.2.3 Put the mattress back on and make sure the bed handle is securely in place.

⚠ **It is important that the bed handle fits tightly against the mattress!**

Now the bed handle "FRIDA" can be used.

- 4.2.4 Mounting the bag

Place the bag in the middle of the handle so that the Velcro strips point inwards and fold them together. Now connect the Velcro strips together.

4.3 Changing the duvet covers or mattress

As a rule, it will not be necessary to remove the "FRIDA" bed handle for this, as it is only on the bed base and has no firm connection to the mattress.

4.4 Instructions for removing the "FRIDA" bed handle (e.g. changing beds or taking it with you when travelling)


Lift or take the mattress down. First, loosen the strap on the end piece and the bed frame. Remove the entire bed handle and place it on a flat surface. Separate the handle by pressing in the safety pins and slowly and evenly pulling the handle out of both base pieces. Detach the base pieces from the end piece by first pressing in the safety pin on one side and slowly pulling out the base piece, then repeat this on the other side.

5.0 Precautions / Warnings

⚠ For safety reasons, it is important that the following conditions are met during assembly and use to avoid the risk of injury or entrapment and, if applicable, the risk of suffocation.

5.1 Dimensions and distances to be observed (based on standard DIN EN 60601-2-52)

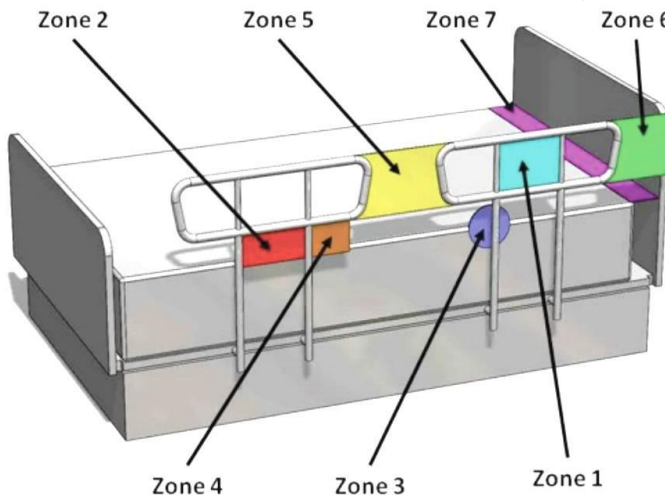Zone 1: / Not applicable to this product.

Zone 2: / Not applicable to this product.

Zone 3: Between the bed handle and the mattress - This area between the inside of the bed handle and the mattress, if too large, can be a risk of head pinching. The maximum allowable distance is 12 cm (4.75 in.).

If possible, however, the bed handle should fit tightly against the mattress.

Zone 4: / Not applicable to this product.

Zone 5: Between Split Bed Handles – If multiple bed handles are used on the same side of the bed, the space between the bed handles can pose a risk of head, neck, or chest pinching. The distance between the bed handles should be at least 32 cm (12.5 in.)

Zone 6: Between the end of the bed handle and the headboard of the bed – The recommended distance between the end of the bed handle and the headboard of the bed is at least 32 cm (12.5 in.), otherwise there is a risk of trapping the head, neck or chest.

Zone 7: Between the headboard or footboard of the bed and the end of the mattress – The size of the mattress should be precisely matched to the bed frame, substructure. If the distance between the headboard or footboard of the bed and the end of the mattress is too large, the risk of trapping increases.

5.2 Mattress features and requirements

There is an increased risk of entrapment and resulting injuries when using mattresses that are too soft or too thin.

To ensure that your mattress is suitable for use with this product, test the firmness of your mattress by pressing a fist-sized object weighing approximately 12 lbs (5.4 kg) onto the edge of the mattress. If the sinking depth is more than 5 cm (2 in.), the mattress used is too soft and the use of the bed rail should be avoided.

The thickness of the mattress should be at least 20 cm (8 in.).

5.3 Warnings, Notices, and Exclusions

- ⚠ The instructions accompanying the product must be followed carefully to ensure that the product is properly installed and constructed. Failure to follow the instructions may result in serious injury to the user(s) of the product.
- ⚠ The product is not intended to support the full weight of a person. This product is only intended as a balance aid and support when sitting down and getting up from bed.
- ⚠ The "FRIDA" bed handle must not be used on cots, bunk beds, waterbeds or inflatable beds.
- ⚠ Under no circumstances should you attach the "FRIDA" bed handle on the handle side to the bed, slatted frame or bed frame. This can lead to serious injury to the user(s) due to entrapment.
- ⚠ Check the positioning of the "FRIDA" bed handle regularly to identify areas where entrapment is possible at an early stage.
- ⚠ With electrically adjustable bed parts, make sure that you do not pinch your limbs when going up and down.
- ⚠ Make sure that children do not play with or hang on to the product.
- ⚠ Additional safety measures should be taken for people who belong to certain risk groups. These include people with a mental impairment (including people who are under the influence of

medication), very tall people or people who are severely overweight (which could change the gap between the FRIDA bed handle and the mattress), as well as people who are generally very restless when awake or asleep.

- ⚠ If you have a medical condition or injury or are taking medications that affect your sense of balance or your ability to stand or sit, please consult your physician before using this product. Improper use of the product may result in serious injury or damage.
- **DISCLAIMER:** Neither the manufacturer nor the distributor is liable for any damage or injury that may result or has resulted from the improper use or assembly of the product.
- In case of installation problems, please contact the selling office (specialist dealer) or directly: rehastage GmbH, Artlandstraße 73, 49610 Quakenbrück / Germany, qm@rehastage.de
- In the event of serious undesirable incidents (deterioration of the user's condition) caused by the use of our product, we are obliged to report it to the Federal Institute for Drugs and Medical Devices. Therefore, we ask you to inform us directly in such a case.

6.0 Cleaning instructions:

Clean the "FRIDA" bed handle regularly with a damp sponge, mild soap and a household disinfectant. Dry it with a dry soft cloth.

Do not use abrasive cleaning agents such as scouring milk or powder, as this can cause destruction or damage to surfaces.

⚠ Warning:

- **Material damage!** Do not use pressure washers, harsh, corrosive chemicals, or abrasive additives for cleaning. After cleaning, make sure that all materials are completely dry to maintain their quality.

Disinfection:

When manually reprocessing a used aid, wipe all parts thoroughly with a surface disinfectant.

Particular care is recommended when disinfecting surfaces with frequent hand/skin contact, such as the handles.

⚠ Warning:

- **Material damage!** When disinfecting, pay attention to the usage and processing instructions of the manufacturer of the cleaning agents and disinfectants.
- The following disinfectants are compatible with the material: Aldehyde-free disinfection based on alcohols (max. 70% isopropyl alcohol)

7.0 Transfer and reuse:

Due to its solid construction, the "FRIDA" bed handle is suitable for reuse. However, the number of reuses depends on the current state of use of the bed handle.

Before being passed on to another user or reused, the FRIDA bed handle must be cleaned, disinfected, inspected for damage and approved by an authorised specialist dealer. This also includes checking the functionally relevant components such as the locks. All technical documentation necessary for safe handling must be handed over.

8.0 Lifespan:

The expected service life of the FRIDA bed handle is usually about 5 years. The prerequisite for this is the intended use, and compliance with the maintenance, cleaning and safety instructions specified in this instruction manual.

However, the service life of the bed handle can be significantly shortened due to disproportionate use, excessive loading, exceeding the maximum weight capacities or lack of care.

9.0 Warranty:

In addition to the statutory warranty, rehastage offers a warranty on this product against manufacturing or material defects of 24 months from the date of purchase. The warranty is only valid against the first user of the item. Keep the invoice that is your proof of purchase.

The warranty is void if the item is repaired or altered by anyone other than rehastage itself or a dealer authorized by rehastage. The warranty obligation expires if repairs or modifications or replacements of individual parts that are not original rehastage spare parts are carried out by a party other than the seller or the supplying company. During the warranty period, rehastage is responsible for the decision to repair the device free of charge or to replace it free of charge.

10.0 Liability:

Rehastage GmbH is only liable if the products are used under the specified conditions and for the specified purposes. We recommend handling the products properly and caring for them according to the instructions. No liability is assumed for damage caused by components and spare parts that have not been approved by Rehastage GmbH. Repairs must only be carried out by authorized dealers or by the manufacturer itself.

11.0 Technical data:

Width:	51 cm
Length:	88 cm
Height of the handle:	adjustable 42 - 53 cm
Net weight:	2.6 kg
Max. user weight:	up to 100 kg
Max. load in use:	Depending on the installation conditions.

12.0 Disposal / recycling of materials:

If the tool was provided to you free of charge, it may not belong to you. If it is no longer needed, please follow the instructions given to you by the organization that provided you with the aid regarding the return. The following section describes the materials used for the aid, with which information on the disposal or recycling of the aid and its packaging. Special regulations regarding disposal or recycling may also apply locally, these must be considered during disposal.

- Steel: Bed Handle "FRIDA"
- Synthetic fibres, plastics: Strap, bag, packaging material
- Carton: Outer packaging

Disposal should be carried out by a waste management company or a public disposal point.

13.0 Explanation of symbols used markings and information to be provided.

	Manufacturer, name and address of the manufacturer.
	Date of manufacture / YYYYMM
	UDI, unique numeric or alphanumeric code for a medical device.
	Authorised representative for medical device manufacturers based outside Switzerland
	Medical device
	Alerts and security measures
	Follow the instructions for use!
	Conformité Européenne: The product complies with the applicable safety, health and environmental requirements of the European Union.

14.0 Accessories Optional:

14.1 Glides for bed handles and adjustable bed side rails

The glides allow easy use of bed rails and bed side rails when using slatted frames with height-adjustable backrests, especially for beds with electrically adjustable slatted frames. Reference: E-glider

14.2 Additional handle "Fridolin"

The additional handle "Fridolin" is mounted in conjunction with the bed handle "FRIDA" and facilitates positioning and straightening in bed as a second handle. Reference: FRIDOLIN 010

The accessories listed here can be purchased from specialist retailers.

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ALBO-HealthCare GmbH